



# Key Topics

- Biochips/Biosensors and Device Fabrication
- Cells, DNA, Proteins
- Micro-fluidics
- Biochip Sensors & Detection Methods
- **Micro-arrays**
- **Lab-on-a-chip Devices**

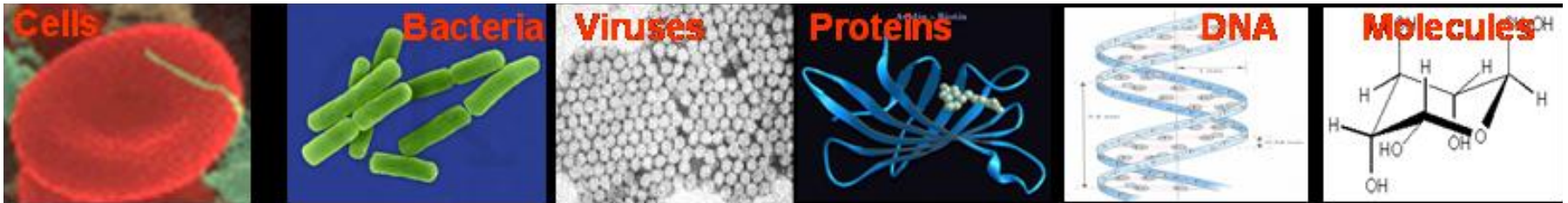

# Optical Detection in Biochips


1. Fluorescence: Markers that emit light at specific wavelengths and enhancement, or reduction (as in Fluorescence Resonance Energy Transfer) in optical signal can indicate a binding reaction
2. Chemiluminescence: Generation of light by the release of energy as a result of a chemical reaction.
  - Light emission from a living organism is termed bioluminescence (sometimes called biological fluorescence),
  - light emission which take place by passage of electrical current is designated electrochemiluminescence.

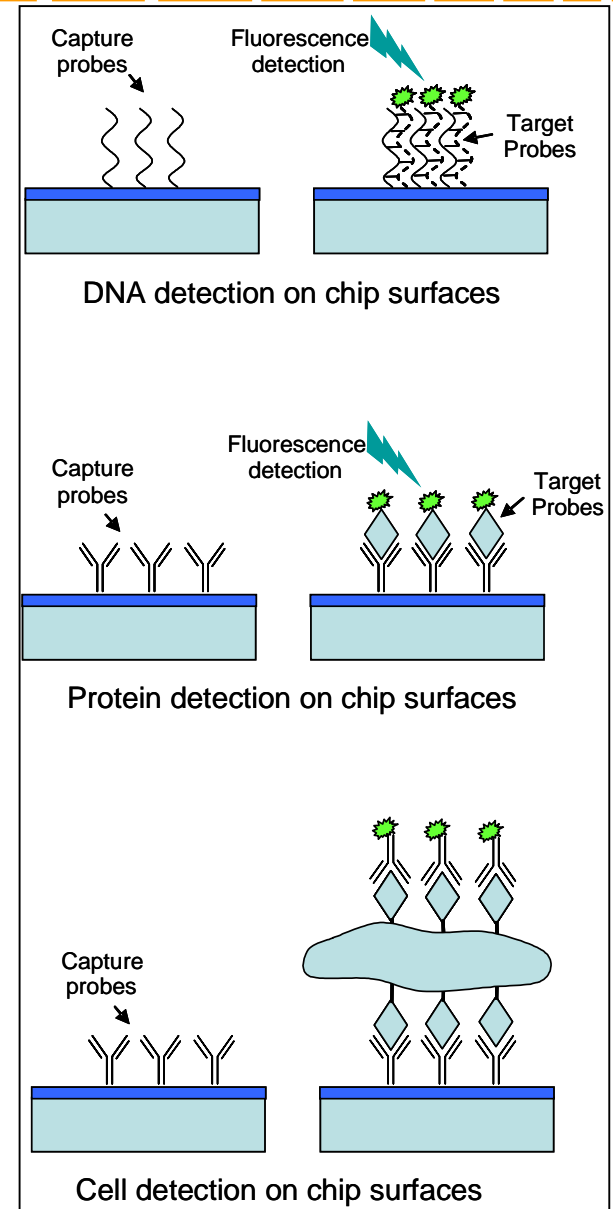


# DNA Hybridization in Microarrays

- Basis for detection of unknown nucleotides
- Example: Bio-chips for identification of DNA
  - Hybridization of an unknown, fluorescently tagged strand with a many known strands - reaction will determine the sequence of the unknown (or vice versa)
  - Strands can be lithographically (Affymetrix) or electronically (nanogen) defined at a specific location

S1	S2	S3
S4	S5	S6
S7	S8	S9

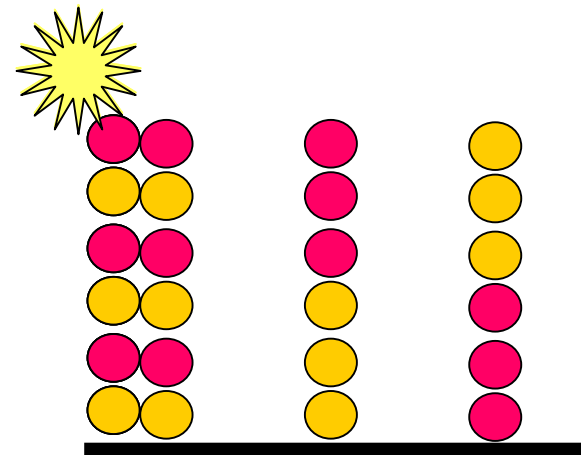

# Electronic Placement of DNA Probes

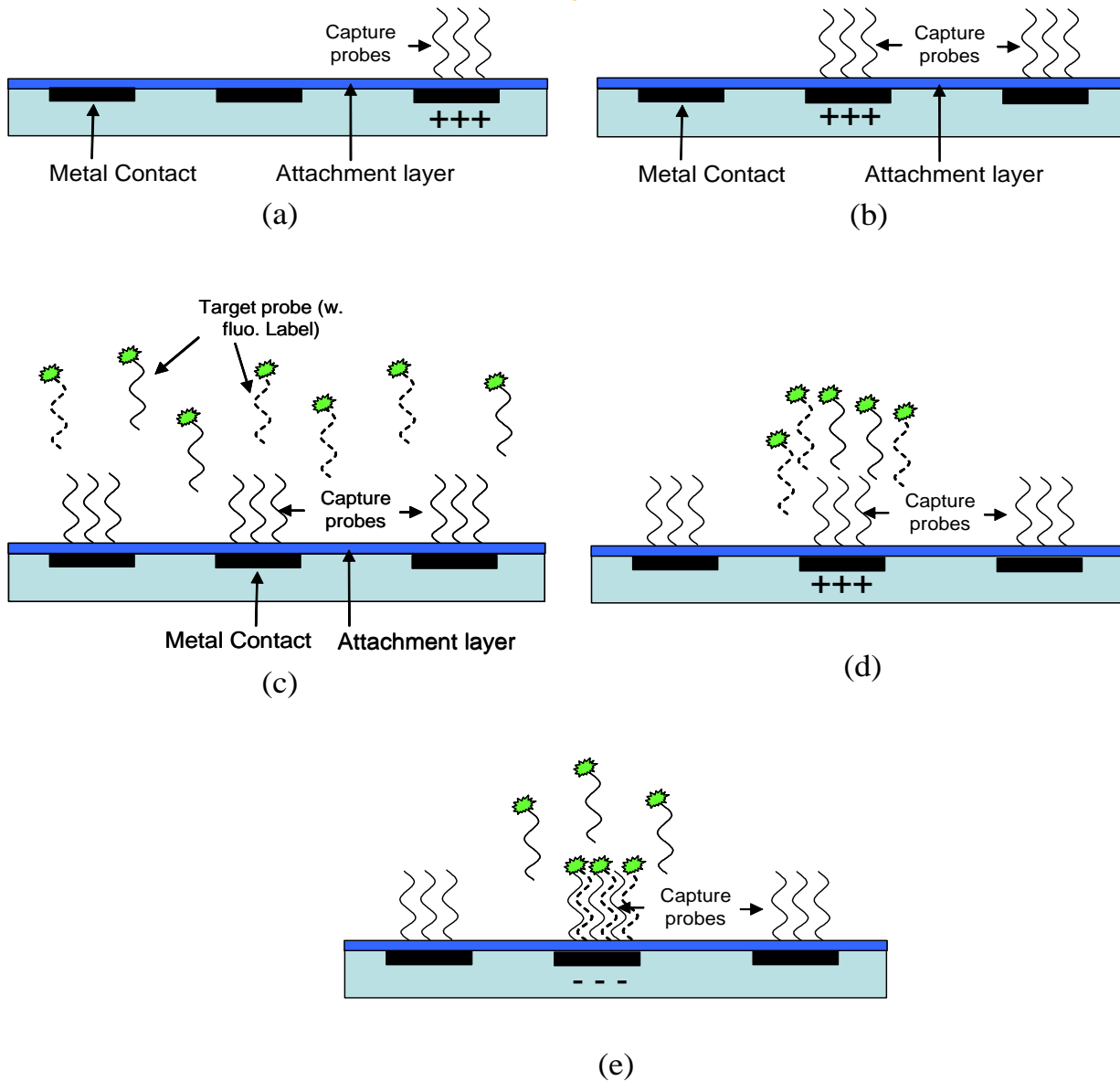

# DNA Biochips (Nanogen)

## Technology Features:




- Biochips for DNA detection, antigen-antibody, enzyme-substrate, cell-receptor and cell separation techniques.
- Takes advantage of charges on biological molecules.
- Small sequences of DNA capture probes to be electronically placed at, or "addressed" to, specific sites on the microchip.




# Technology Features

## Hybridization.


- A test sample can be analyzed for the presence of target DNA molecules by determining which of the DNA capture probes on the array bind, or hybridize, with complementary DNA in the test sample.
- **Fluorescence output**

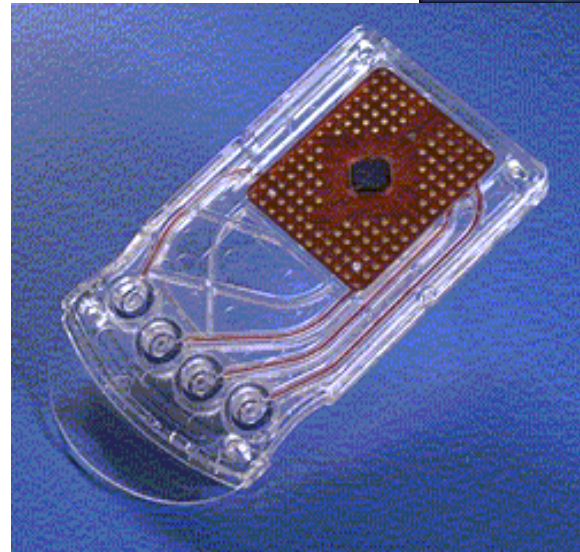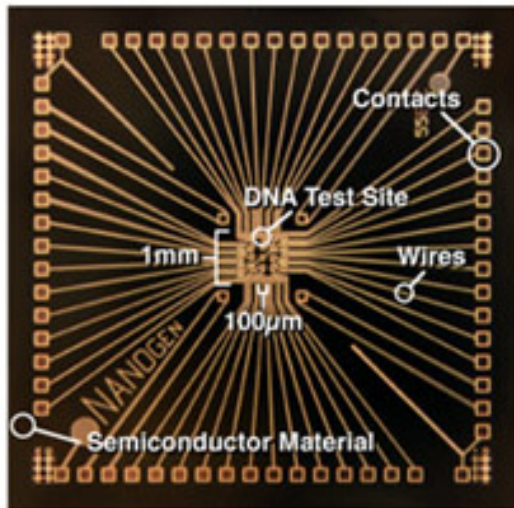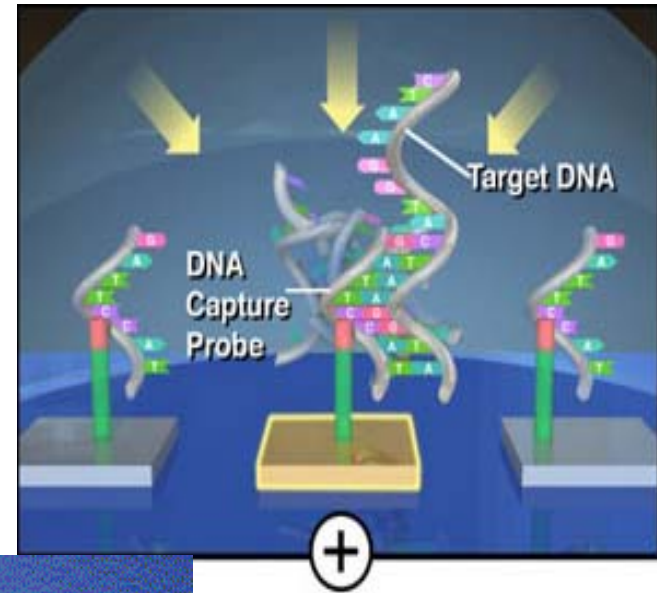

# Light Directed DNA Synthesis on a chip (Affymetrix)

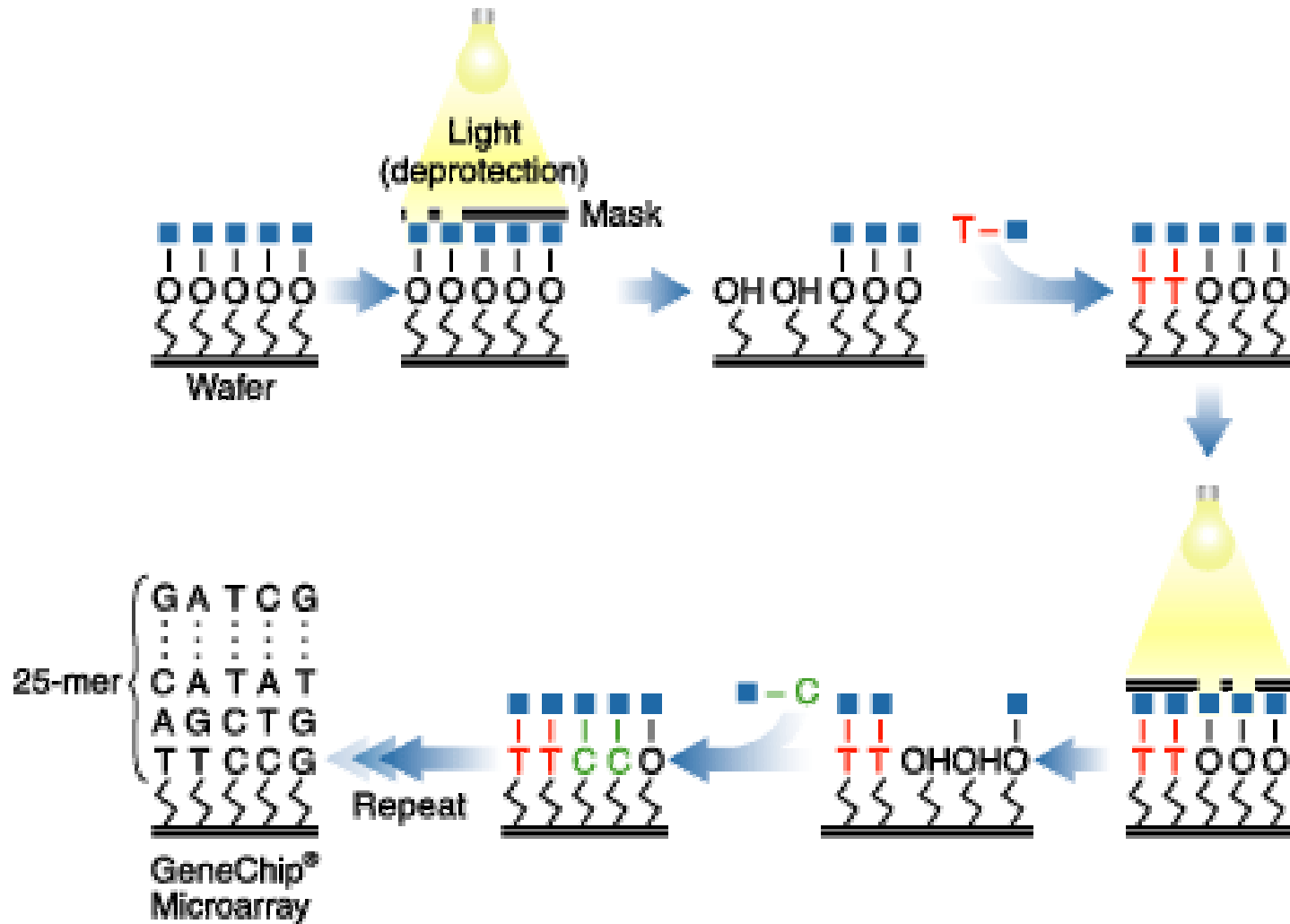

# Light Directed DNA Synthesis on a chip (Affymetrix)

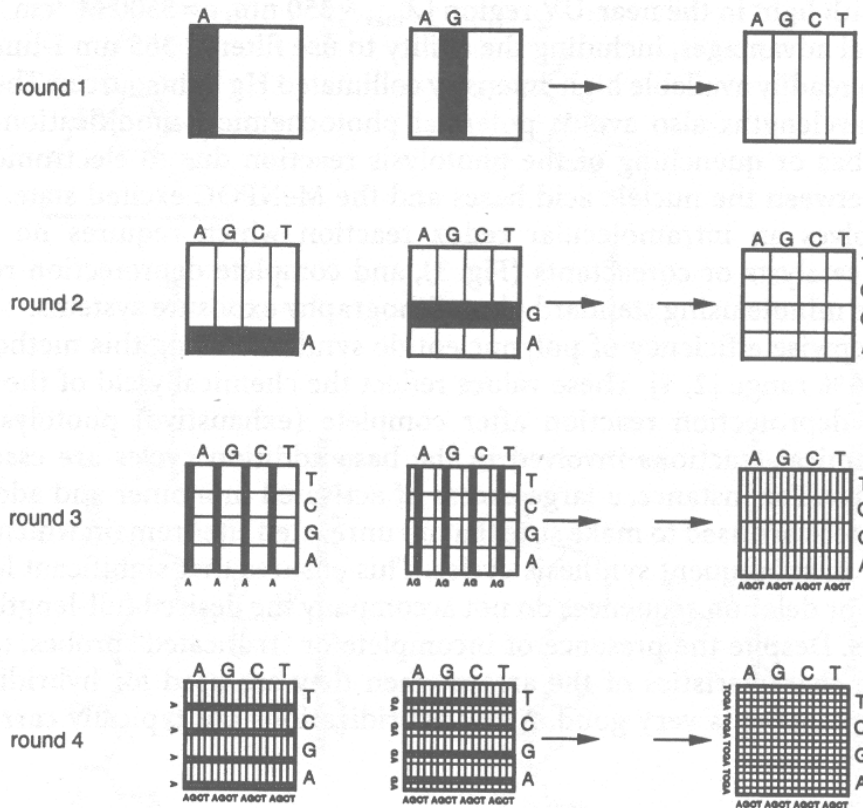


**Table 1.** Combinatorial synthesis of polynucleotide probe arrays


Probe Length	Chemical Steps	Number of Possible Probes
4	16	256
8	32	65 536
12	48	16 777 216
16	64	$\sim 4.3 \times 10^9$
20	80	$\sim 1.1 \times 10^{12}$



# Light Directed DNA Synthesis on a chip (Affymetrix)


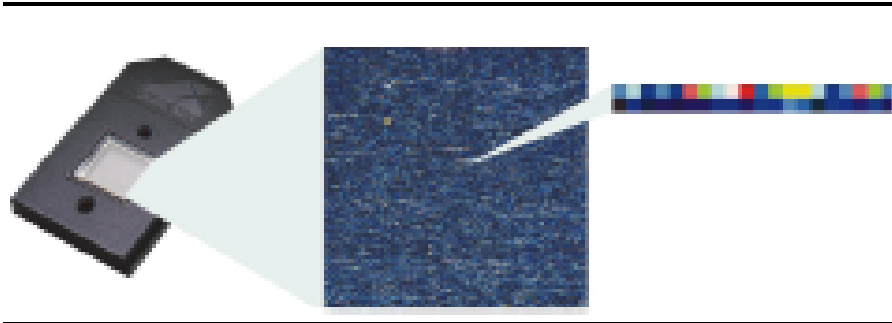


- Fluorescence detection
- Ultimately will limit size of pixel in array

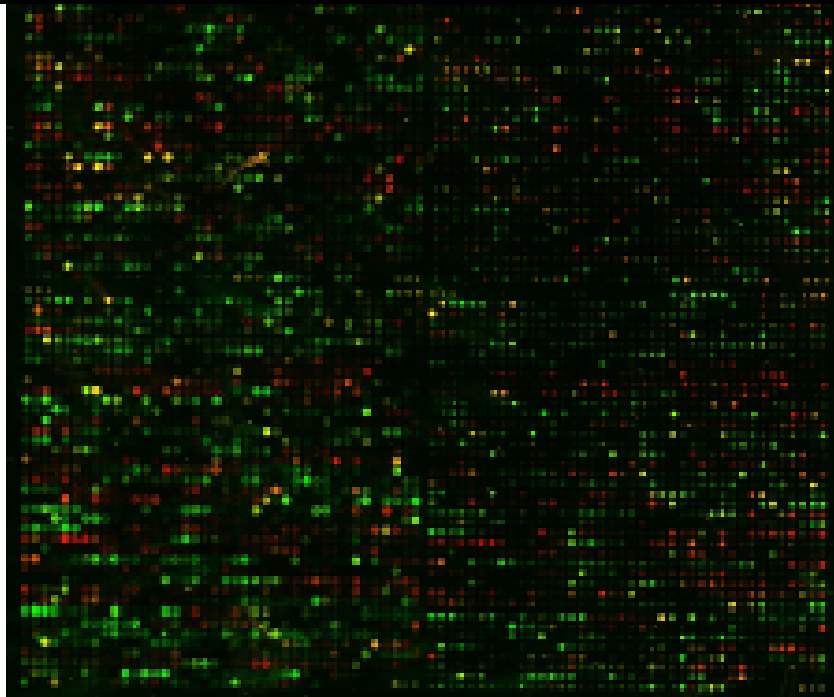

**Applications:**  
Polynucleotide array  
HIV resequencing  
mRNA expression monitoring

# Protein Arrays

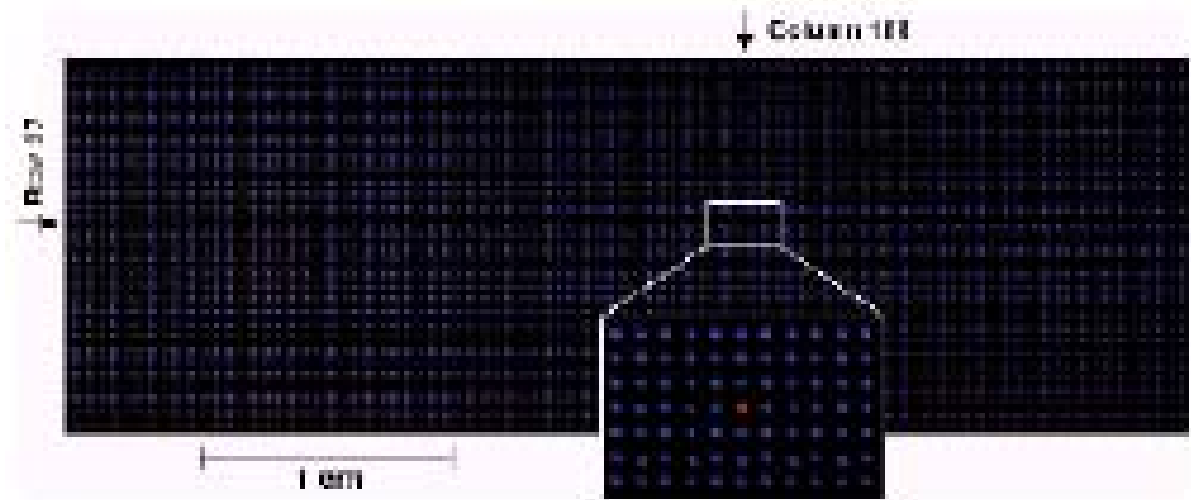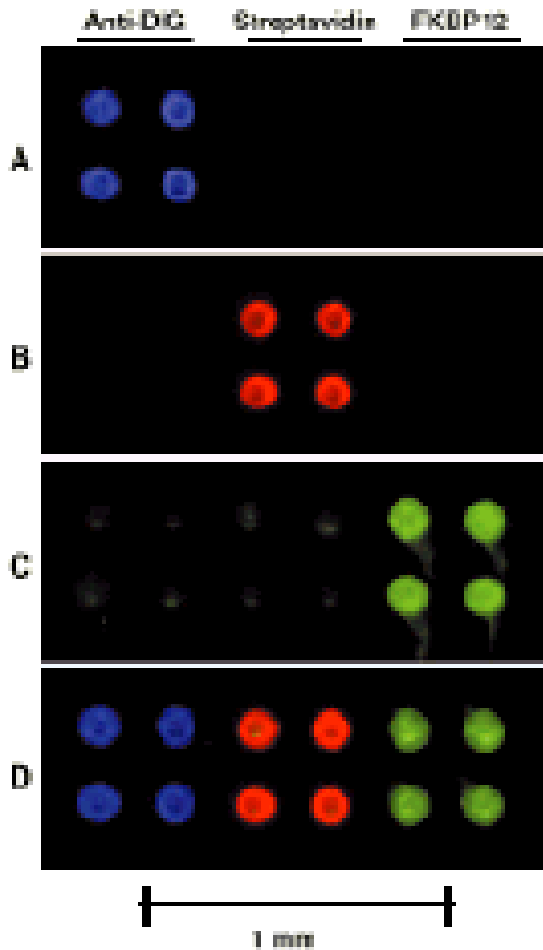

- Protein-Protein Interactions
- Protein small molecule interactions
- Derivatized substrates – glass, plastics
- High Throughput screening of chemical compounds



# Note: Sensor Arrays

---

- Any of the individual sensors described earlier can be used in an array format to make micro/nano sensor arrays.
- The sensors in the array need addressing
- Each sensor can be functionalized with different bio-receptor molecule to detect different entities
- Examples, cantilever array, electrochemical detection in electrode arrays, cellular arrays for chemical detection, etc.


# Lab-on-a-Chip/Integrated Devices

- Single chip device for DNA electrophoresis
- Sample loading and metering
- PCR on a chip (faster temperature cycling due to reduced thermal mass)
- Gel electrophoresis on chip

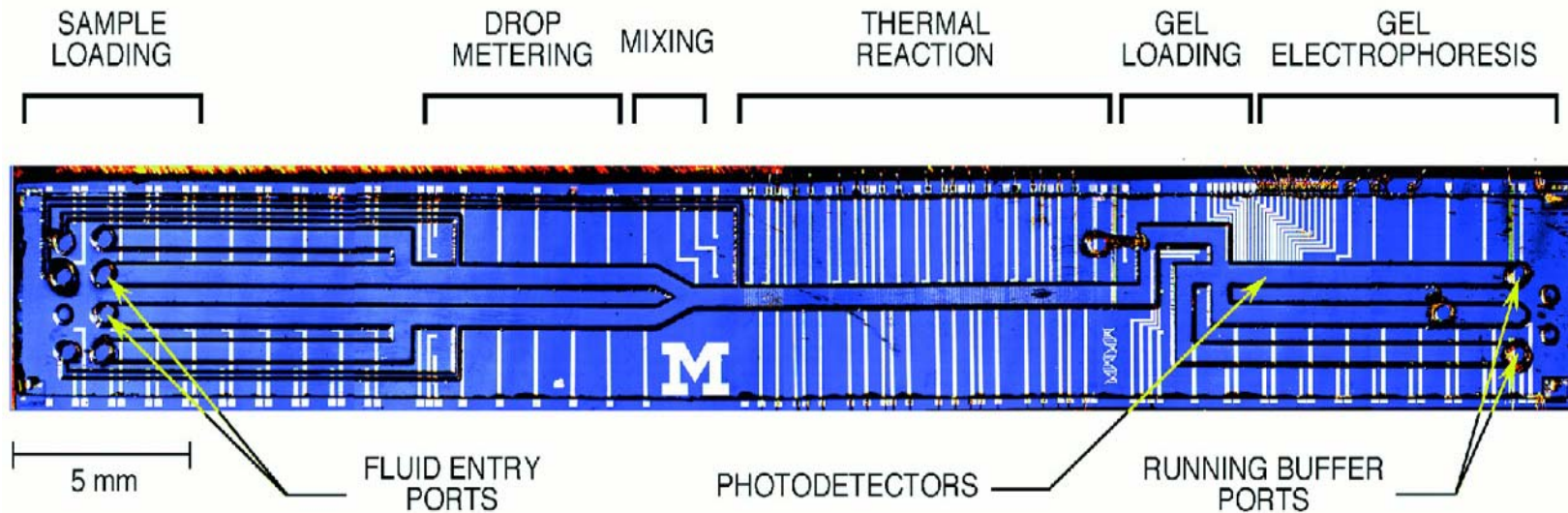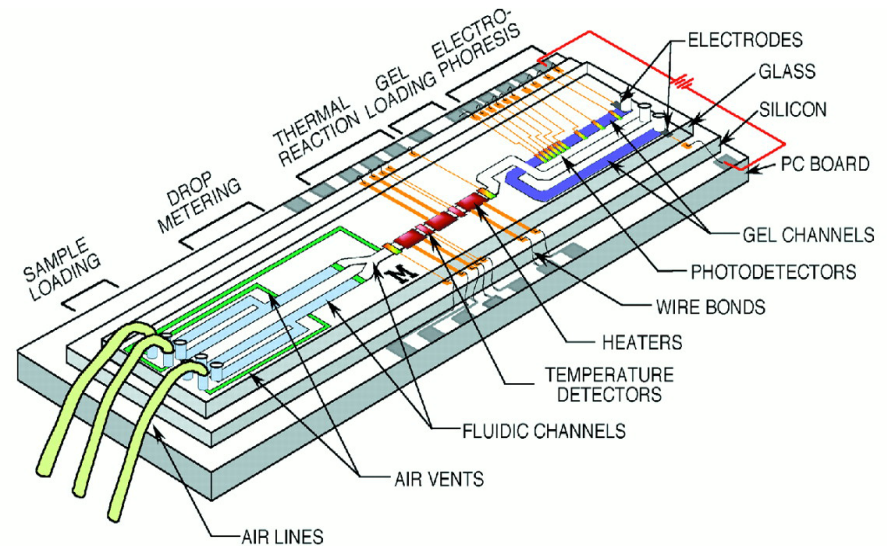

# CD Format Biochips



- Micro-fluidic devices on a CD type platform using centrifugal and capillary forces for liquid transport
- Cheap plastic CDs
- Optical detection systems

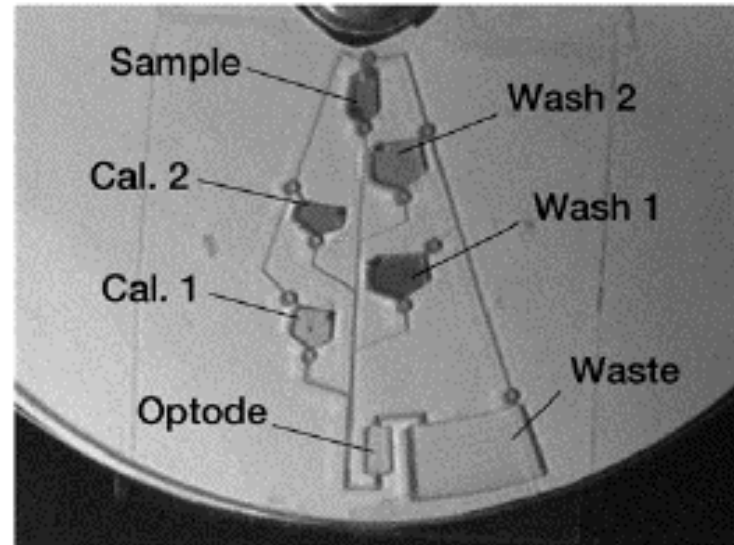


Flow order: Cal. 1 → Wash 1 → Cal. 2 → Wash 2 → Sample

# Cellular Analysis on Chip



- Plastic biochips using hydrodynamic transport of cells
- Electric field mediated lysing
- Fluorescence detection (off-chip detectors)
- Analysis time of about 10 cells/minute

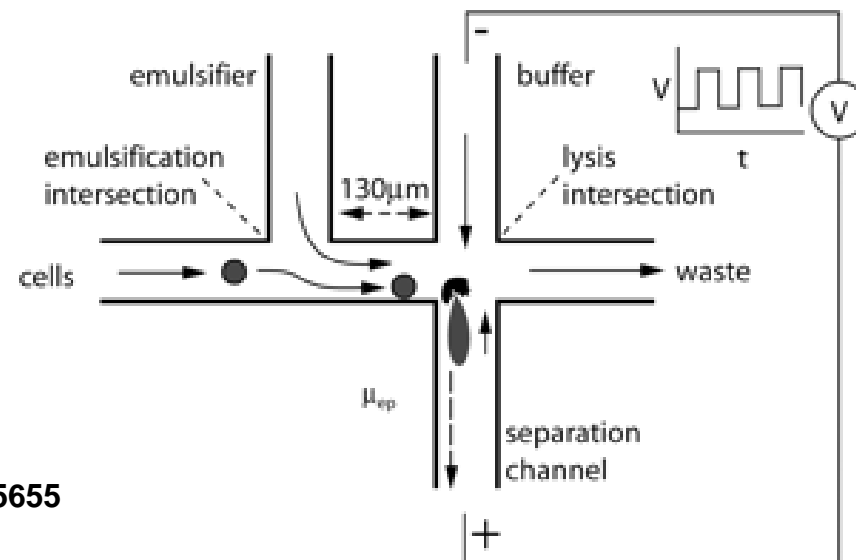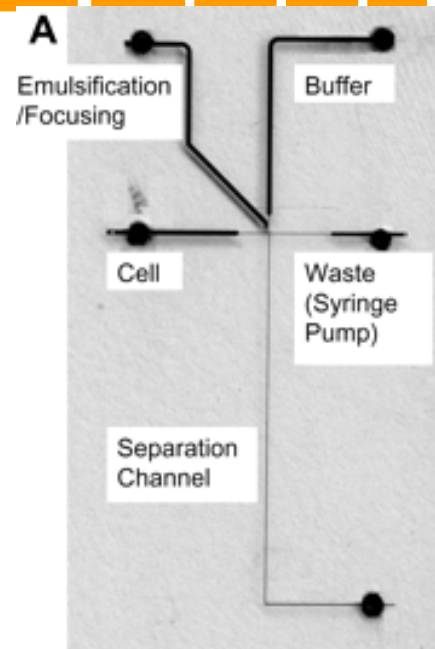


# Polymer $\mu$ Sensor and Actuator

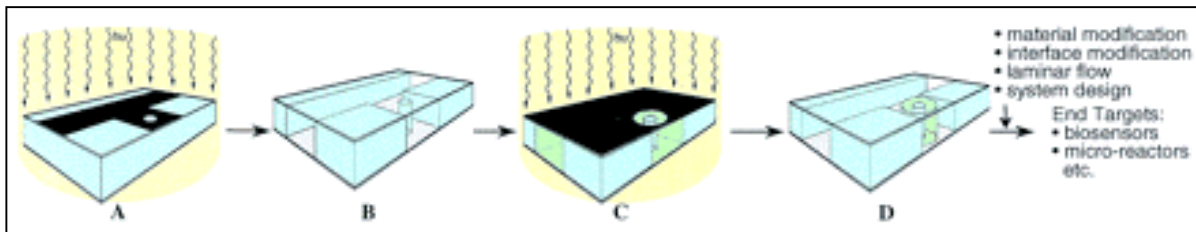

Process flow for the preparation of a hydrogel valve.



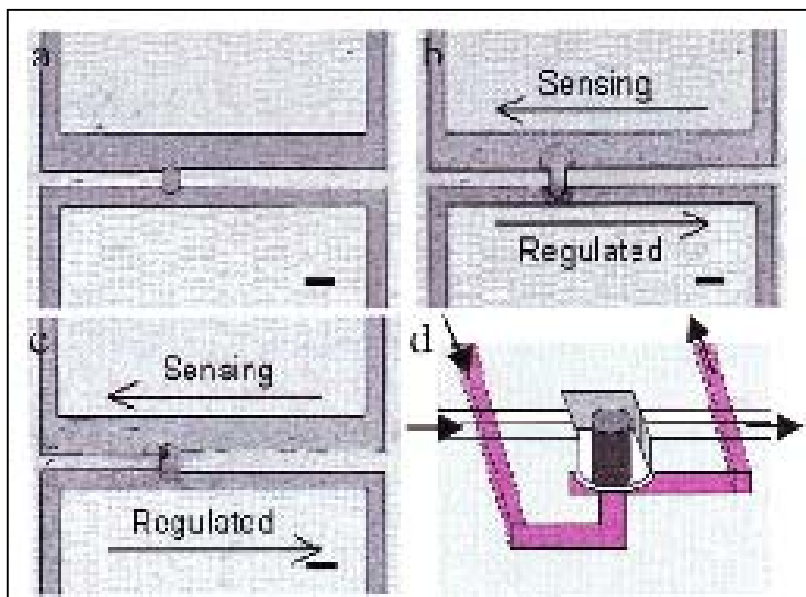
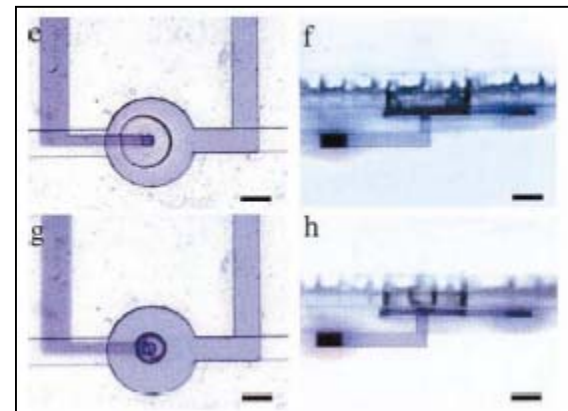


Hydrogel valve designs (2D and 3D)

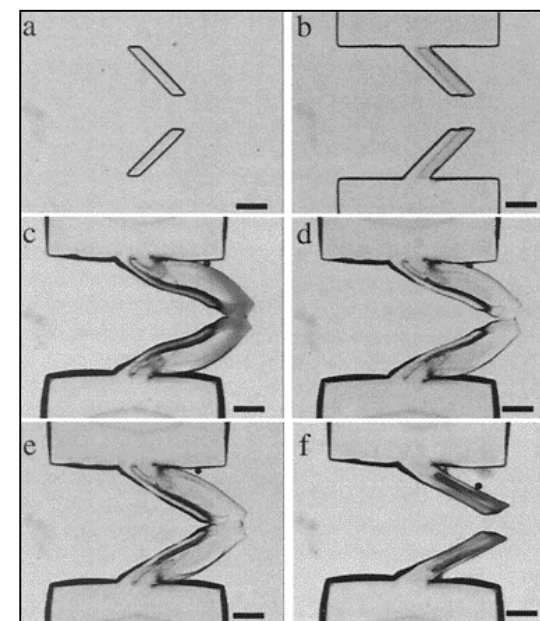

A biomimetic valve based on bistrif hydrogel.


# DNA Capillary Electrophoresis

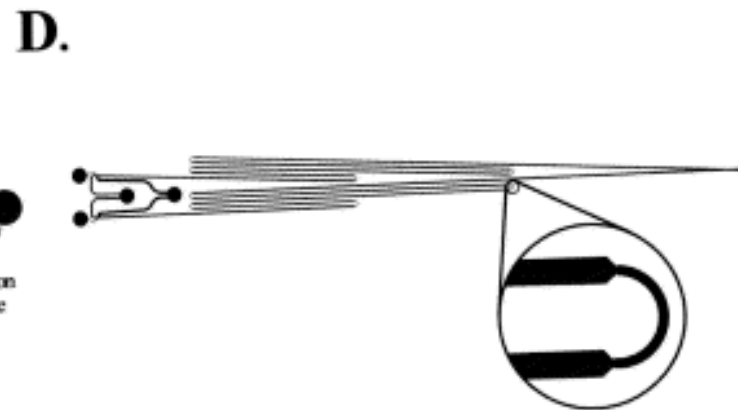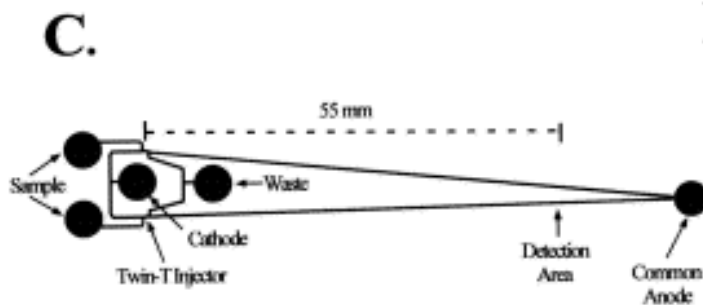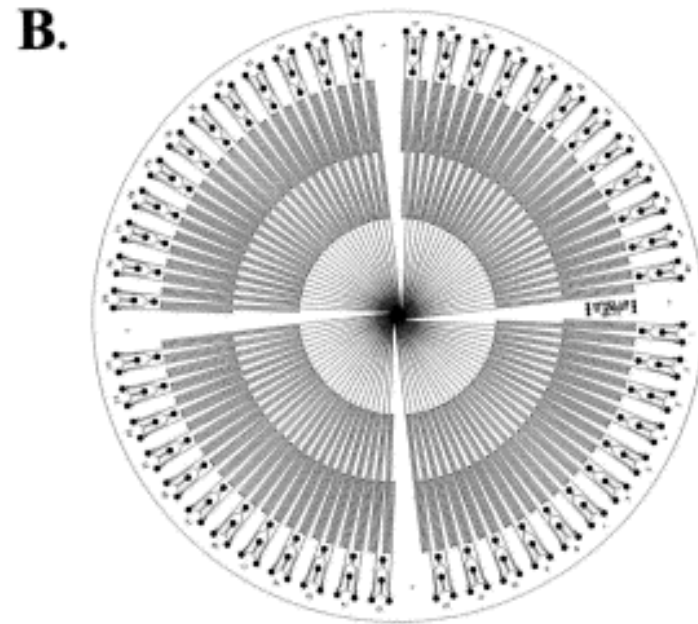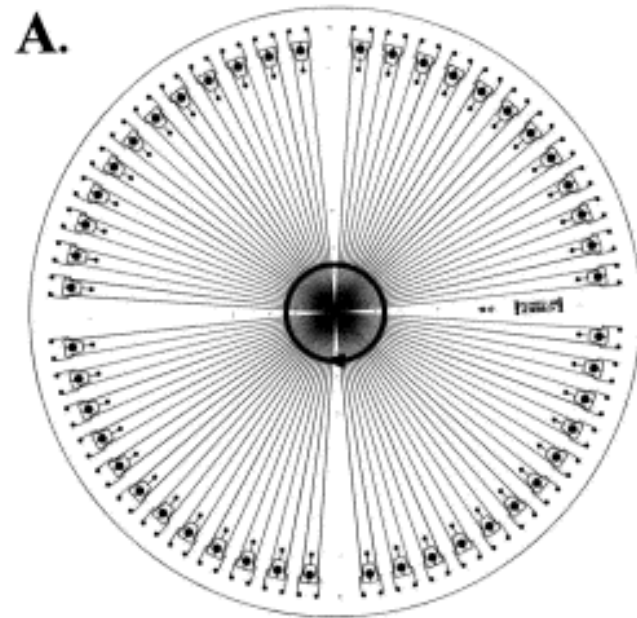

**Design of the 96-channel CAE microplate and radial scanner.  
Mask pattern used to form the 96 straight channel  
radial microplate on a 150-mm diameter wafer.**






# Integrated Systems for Study of Microorganisms and Cells

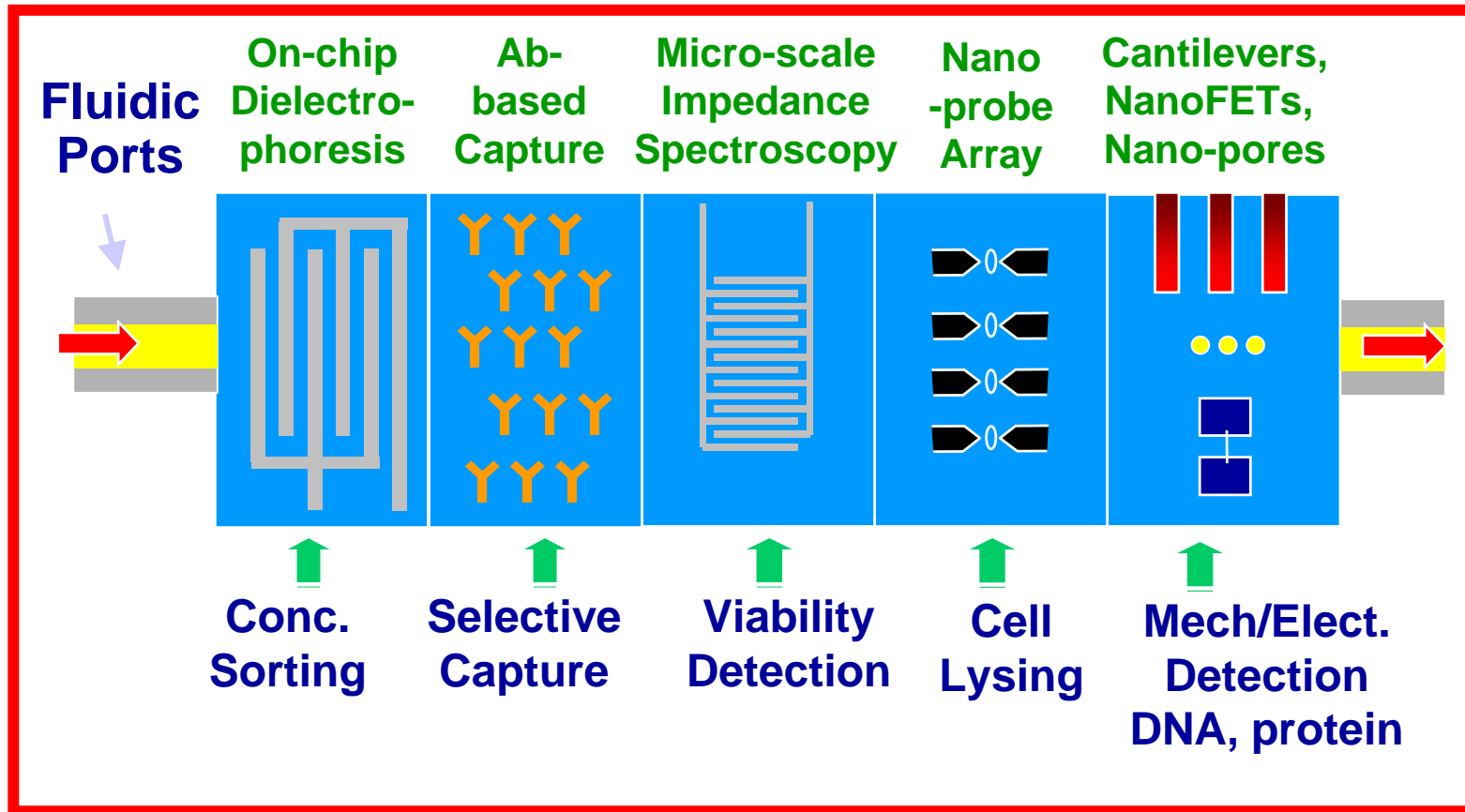

“Lab on a Chip” for  
Enabled by BioMEMS and  
Bionanotechnology

# Micro-fluidic Polymer Devices for Culture Bacteria and Spores


- Growth of bacteria inside a micro-fluidic polymer chip
- Rapid detection and reduced time to result

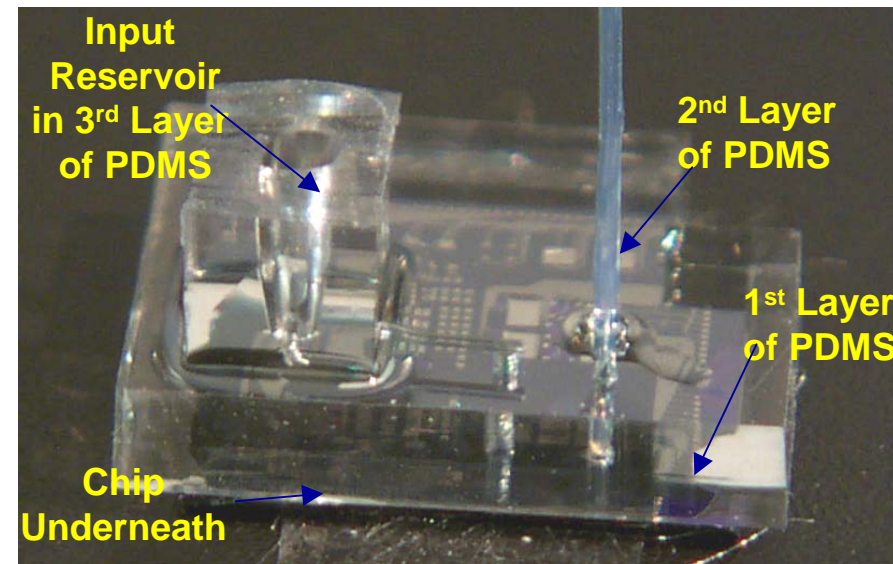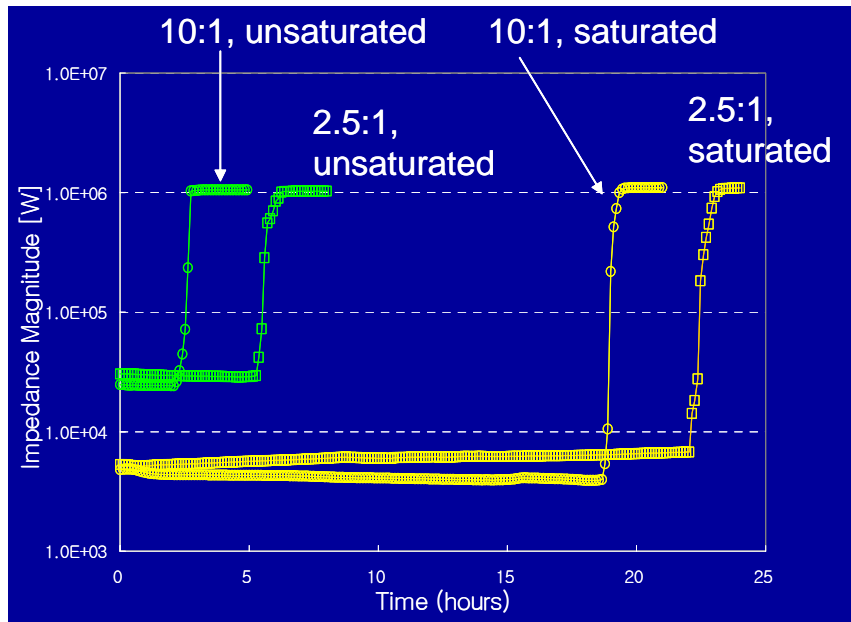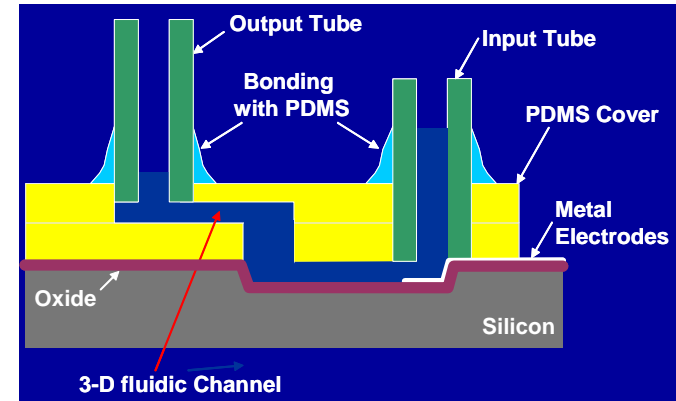


Silicon Base, 3 PDMS layers, Top I/O port

# Future Directions

- Integrated device for analysis of single cells – applications and fundamental science
- Building cell by cell/tissue engineering using micro and nano fabrication techniques
- Integrated diagnostics and therapeutics (drug delivery)
- Tools for genetic manipulation of microorganisms and viruses – synthetic biology

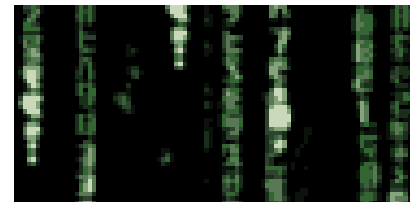

# Acknowledgements

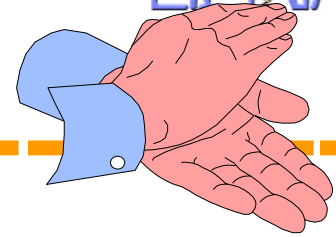

## Research Scientists/Post-docs:

- Dr. Demir Akin
- Dr. Dallas Morissette
- Dr. Rafael Gomez

## Graduate Students:

- Sangwoo Lee
- Haibo Li
- Amit Gupta
- Hung Chang
- Yi-Shao Liu
- Samir Iqbal
- Oguz Elibol
- Angelica Davilia
- Kidong Park

## Industries:

- BioVitesse, Inc. Co-Founder

## Funding Agencies

- US Department of Agriculture (Food Safety Engineering Center)
- NASA Institute on Nano-electronics and Computing
- NSF, NSF Career Award
- National Institute of Health
- DARPA Nanotechnology Research
- Discovery Park at Purdue University

## Faculty Collaborators

- Prof. D. Bergstrom (Med Chem)
- Prof. A. Bhunia (Food Science)
- Prof. M. Ladisch (Ag& Bio Engr)

## Special Thanks

- Prof. S. Broyles (BioChem)
- Profs. D. Datta, D. Janes (ECE, NASA INAC), J. Cooper (BNC)

